

Pro-Tip[®] *TURBO*



bringing
air and water
to life

YOUNG[®]
INNOVATIONS

designed to keep

air and water

separate

Many dental procedures require instant, clean, dry air. Pro-Tip Turbo has a unique 2-chamber system to separate air from water. This guarantees instant dry air when and where you need it. Pro-Tip Turbo also produces a perfect mist every time.



- High power air jet
- Prevents cross contamination
- Dry air every time and perfect misted spray
- Locking system secures tip safely
- Sealed separation of air and water
- Rigid for effective cheek retraction
- Wide range of converters to fit most syringes

"I have used disposable 3-in-1 tips for a number of years now without being fully satisfied. However, clean, dry air is delivered predictably everytime with Pro-Tip Turbo - essential for today's adhesive restorative dentistry. They certainly live up to the 'turbo' name with air water being delivered at speed. The adapter was easy to fit to my A-dec unit and the bayonet design prevents any tips shooting off. They offer great value too always a consideration for disposable items."







Dr. James Robson BDS

how it works

Pro-Tip[®] TURBO

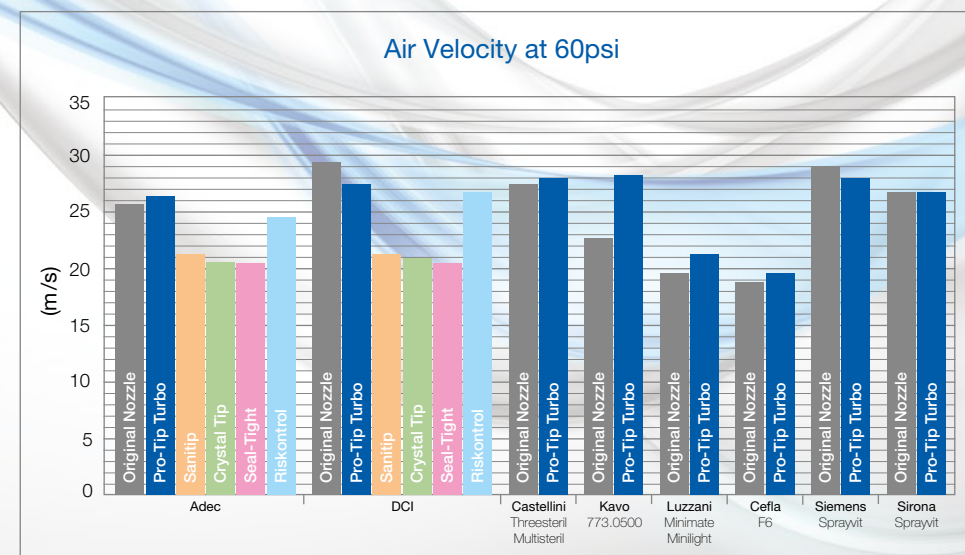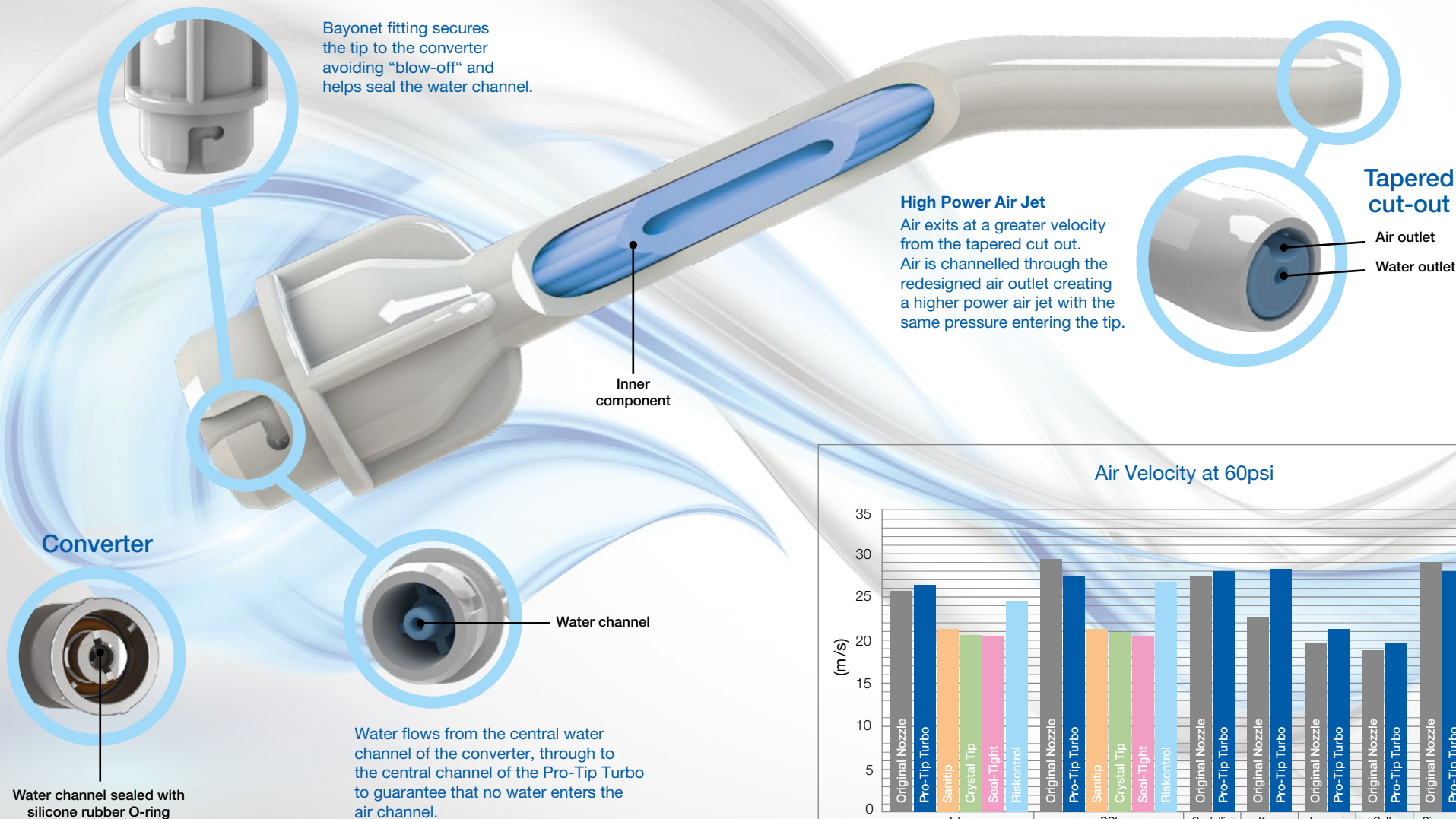


Choose a converter to fit your syringe...

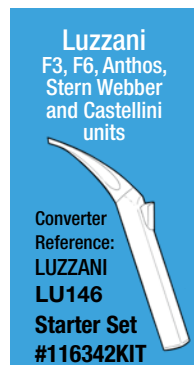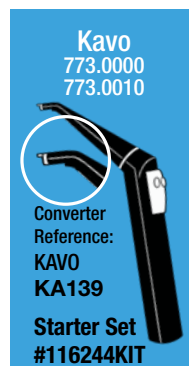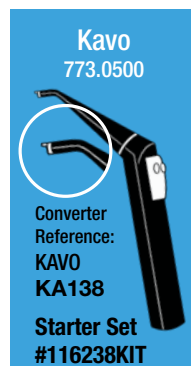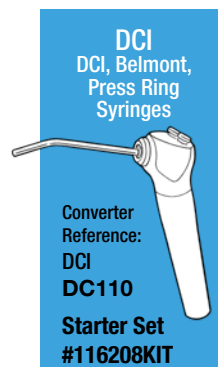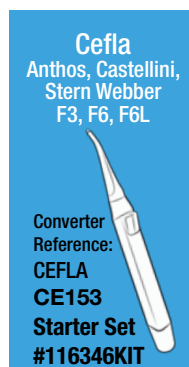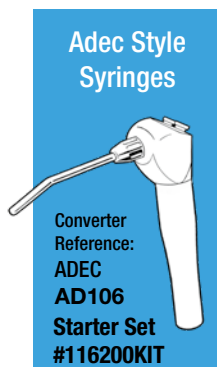

Easy installation

1. Remove existing syringe nozzle and/or converter then attach Pro-Tip Converter to syringe.
2. Push Pro-Tip onto converter. Twist Pro-Tip whilst holding the outer sleeve of the converter to engage locking feature.
3. Pro-Tip is securely attached to converter when the ribs are aligned.

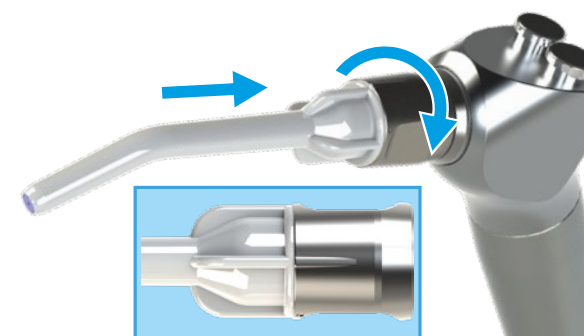

Visit our website

www.protipdental.eu

to view the full range of converters
and installation videos.

YOUNG[®]

INNOVATIONS



Young Innovations Europe GmbH

Mittermaierstr. 31 | D-69115 Heidelberg | Germany

IE: +353 (58) 45 966 | info@ydnt.eu | www.ydnt.eu

For more information visit: www.protipdental.eu